



Washing and
disinfection in the
hospital sector

2018

Machine accessories

Water Treatment: Demineralizer

WO-10 (art. 090008)

PURIFIER FOR REVERSE OSMOSIS IN STAINLESS STEEL

- electronic control panel
- alarm for filter depletion
- alarm for lack of water supply
- alarm for anti-flood with automatic blocking system
- automatic timed wash of membranes in case of prolonged inactivity
- counters for hours and liters
- water outlet conductivity (max): 5-30 $\mu\text{S}/\text{cm}$
- hourly production: 120 liters/h
- production and conductivity depend on both the chemical characteristics of the water inlet and the clogging of the membranes
- recovery rate: 40%
- pre-treatment on water inlet: activated carbon filter
- water inlet pressure (min / max): 1,5 / 6 bar
- water inlet hardness (max): 40 °f
- water inlet temperature (max): 30 °C
- water inlet conductivity (max): 2000 $\mu\text{S}/\text{cm}$
- external dimensions LxDxH: 143 x 450 x 456 mm
- weight (net / gross): 15 / 17 kg
- electrical connection: 1/N/PE 230 V - 50 Hz - 150 W



WP3000 (art. 900187)

WATER PURIFIER - MIXED BED RESIN

- to provide demineralized water from tap water
- conductivity sensor included to indicate exhausted resins
- external dimensions LxDxH: 300 x 660 x 850 mm
- complete with CL3000 column
- resins kit MI700WP
- auxiliary draw-off valve



Water Treatment: Softener

WS11EB (art. 903196)

PROFESSIONAL WATER SOFTENER - RESIN CAPACITY 11 LITRES

- timed 5 phase regeneration
- hardness regulation mixer
- cycling capacity of water production: 1650 liters at 0 °f with incoming water at 40 °f
- salt consumption for regeneration cycle 1,0 kg about
- plastic reinforced by incorporated fiberglass
- dimensions LxDxH: 300 x 470 x 660 mm
- filter pot holder of 9" included
- 3/4" IN and OUT connections



WS14EB (art. 903197)

PROFESSIONAL WATER SOFTENER - RESIN CAPACITY 14 LITRES

- timed 5 phase regeneration
- hardness regulation mixer
- cycling capacity of water production: 2250 liters at 0 °f with incoming water at 40 °f
- salt consumption for regeneration cycle 1,2 kg about
- plastic reinforced by incorporated fiberglass
- dimensions LxDxH: 300 x 470 x 815 mm
- filter pot holder of 9" included
- 3/4" IN and OUT connections



WS17EB (art. 903198)

PROFESSIONAL WATER SOFTENER - RESIN CAPACITY 17 LITRES

- timed 5 phase regeneration
- hardness regulation mixer
- cycling capacity of water production: 2700 liters at 0 °f with incoming water at 40 °f
- salt consumption for regeneration cycle 1,8 kg about
- plastic reinforced by incorporated fiberglass
- dimensions LxDxH: 300 x 470 x 1070 mm
- filter pot holder of 9" included
- 3/4" IN and OUT connections



Trolleys

Base trolleys

CBS1 (art. 900785)

BASIC TROLLEY FOR 2ND/3RD/4TH LEVELS

- sprayer washing
- incorporated sprayer
- for holding baskets and supports
- made of stainless steel



CBS2 (art. 900787)

BASIC TROLLEY FOR 1ST LEVEL

- sprayer washing
- bottom tank sprayer
- for holding baskets and supports
- made of stainless steel



CS1-1 (art. 901696)

UPPER TROLLEY WITH SPRAYER

- suitable for holding of supports and baskets
- usable area LxD: 460 x 450 mm
- usable height adjustable on two positions: 190 mm or 240 mm
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

CS2 (art. 900032)

LOWER TROLLEY

- suitable for holding of supports and baskets
- usable area LxD: 480 x 500 mm
- usable height: 230 mm
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

CS3 (art. 900323)

UNIVERSAL SUPPORT TROLLEY

- n. 2 washing levels
- upper level with incorporated sprayer arm
- lower level with machine bottom sprayer
- capacity 4 baskets, codes CB18 - CB18/100 - CB18/120 – CPS20
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

CSKSTOV (art. 902815)

UPPER AND LOWER TROLLEY FOR THE TREATMENT OF 14 COVERS

- suitable for washing and thermal disinfection of table covers
- made of anti corrosion coated stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

D-CS1 (art. 901204)

UPPER TROLLEY WITH SPRAYER 45CM

- suitable for holding of supports and baskets
- usable area LxD: 310 x 430 mm
- usable height: 250 mm
- made of stainless steel



Machines

WD2145

D-CS2 (art. 901205)

LOWER TROLLEY 45 CM

- suitable for holding of supports and baskets
- usable area LxD: 330 x 450 mm
- usable height: 250 mm
- made of stainless steel



Machines

WD2145

MAGBAR10 (art. 903457)

MAGNET GROUP FOR AUTOMATIC TROLLEY RECOGNITION

- to be installed on the washing trolley
- allows the automatic trolley recognition
- to be used with RACKREC60 optional

Machines


WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM	WD7010	WD7010M	WD7010T	WD7010TM

Handpieces and cannules trolleys

KITD-SM4E (art. 903961)

WASHING BASKET FOR HANDPIECES

- CSKDIN basket with collector for injection washing of turbines
- supplied with connection set composed of 4 SM1E supports

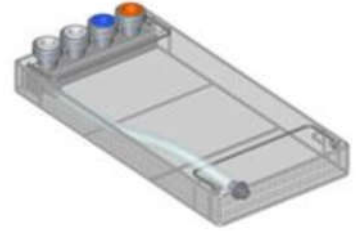

Surgical instruments trolleys

CSK1010 (art. 903368)

5 LEVELS WASHING RACK

- suitable for washing specific supports and/or baskets holding instruments
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- capacity: 10 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM	WD7010	WD7010M	WD7010T	WD7010TM

CSK1010FLEX (art. 903662)

5 LEVELS FLEXIBLE WASHING RACK

- multipurpose rack for simultaneous washing of standard and MIC surgical instruments
- suitable for washing specific supports and/or baskets holding instruments KITMIC10 (max. 4)
- 2, 4, 5, removable levels for washing different instrument heights
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- n. 4 connections for injection accessories (not included)
- capacity: 10 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM	WD7010	WD7010M	WD7010T	WD7010TM

CSK1210 (art. 903369)

6 LEVELS WASHING RACK

- suitable for washing specific supports and/or baskets holding instruments
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- capacity: 12 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM	WD7010	WD7010M	WD7010T	WD7010TM

CSK1210FLEX (art. 903408)

6 LEVELS FLEXIBLE WASHING RACK

- multipurpose rack for simultaneous washing of standard and MIC surgical instruments
- suitable for washing specific supports and/or baskets holding instruments KITMIC10 (max. 3)
- 2, 3, 5, 6 removable levels for washing different instrument heights
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- n. 3 connections for injection accessories (not included)
- capacity: 12 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM	WD7010	WD7010M	WD7010T	WD7010TM

CSK1210UNIV (art. 903453)

FLEXIBLE WASHING RACK

- multipurpose rack for washing of standard instruments, MIC instruments, OP Shoes and containers
- suitable for washing specific supports and/or baskets holding instruments KITMIC10 (max. 3)
- spraying arms on chamber ceiling and bottom for upper and lower levels
- provided with intermediate level with integrated spraying arm and upper level
- n. 3 connections for injection accessories (not included)
- additional PLAV10 levels with integrated spraying arm are not included (max. 3)
- supports and accessories (not included) can be additionally installed
- made of stainless steel



Machines

WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM	WD7010	WD7010M	WD7010T	WD7010TM

CSK1515 (art. 903401)

5 LEVELS WASHING RACK

- suitable for washing specific supports and/or baskets holding instruments
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- capacity: 15 CSKDIN baskets (not included)
- made of stainless steel

CSK1515FLEX (art. 904181)

5 LEVELS FLEXIBLE WASHING RACK

- multipurpose rack for simultaneous washing of standard and MIC surgical instruments
- suitable for washing specific supports and/or baskets holding instruments KITMIC10 (max. 4)
- 2, 4, 5, removable levels for washing different instrument heights
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- n. 4 connections for injection accessories (not included)
- capacity: 15 CSKDIN baskets (not included)
- made of stainless steel

CSK1815 (art. 904071)

6 LEVELS WASHING RACK

- suitable for washing specific supports and/or baskets holding instruments
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- capacity: 18 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD7015

WD7015T

CSK1815FLEX (art. 903863)

6 LEVELS FLEXIBLE WASHING RACK

- multipurpose rack for simultaneous washing of standard and MIC surgical instruments
- suitable for washing specific supports and/or baskets holding instruments KITMIC10 (max. 4)
- 2, 3, 5, 6 removable levels for washing different instrument heights
- spraying arms on chamber ceiling and bottom for upper and lower levels
- intermediate levels with integrated spraying arm
- n. 4 connections for injection accessories (not included)
- capacity: 15 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD7015

WD7015T

CSK4B (art. 902802)

2-LEVEL WASHING TROLLEY

- suitable for decontamination of specific support and/or baskets holding instruments
- upper level with spray arm
- lower level directly above the tank spray arm
- capacity: 4 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

CSK4BMIC (art. 903750)

2-LEVEL WASHING TROLLEY

- suitable for decontamination of specific supports and/or baskets holding instruments
- suitable for cleaning cannulated instruments with basket KITMIC10 (not included)
- upper level with spray arm
- lower level directly above the tank spray arm
- capacity: 4 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD1160
WD4190

WD1160P4
WD4190P4

WD3060
WD5090T

WD4060
WD5090TFO

WD4060P4

CSK6B (art. 902999)

3-LEVEL WASHING TROLLEY

- suitable for decontamination of specific support and/or baskets holding instruments
- upper levels with spray arm
- lower level directly above the tank spray arm
- capacity: 6 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD1160
WD4190

WD1160P4
WD4190P4

WD3060
WD5090T

WD4060
WD5090TFO

WD4060P4

CSK8B (art. 902803)

4-LEVEL WASHING TROLLEY

- suitable for decontamination of specific support and/or baskets holding instruments
- upper level with spray arm
- lower level directly above the tank spray arm
- capacity: 8 CSKDIN baskets (not included)
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4


WD5090T

WD5090TFO

PLAV10 (art. 903405)

SINGLE WASHING LEVEL

- intermediate level with integrated spraying arm
- capacity: 2 CSKDIN baskets (not included)
- made of stainless steel

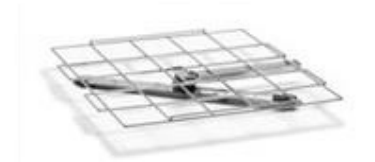

Machines

WD6010

WD6010E

WD6010EM

WD6010M

WD6010T

WD6010TM

WD7010

WD7010M

WD7010T

WD7010TM

WD-TRST (art. 902467)

MAGNET GROUP FOR AUTOMATIC TROLLEY RECOGNITION

- to be installed on the washing trolley
- allows the automatic trolley recognition
- to be used with WD-TRS

Operating room shoes trolleys

D-ZCS1 (art. 901894)

UPPER LEVEL TROLLEY FOR OPERATING ROOM SHOE WASHING (45CM MACHINES)

- equipped with sprayer arm
- capacity: 3 pairs of standard operating room shoes
- made of stainless steel



Machines

WD2145

D-ZCS2 (art. 901895)

LOWER LEVEL TROLLEY FOR OPERATING ROOM SHOE WASHING (45CM MACHINES)

- level with machine bottom sprayer
- capacity: 4 pairs of standard operating room shoes
- made of stainless steel



Machines

WD2145

SZ1010 (art. 903370)

SHOE WASHING SUPPORT STRUCTURE

- capacity: 10 pairs of standard OP shoes
- can be placed on rack CSK1210UNIV (max. 2)
- made of stainless steel



Machines

WD6010

WD6010E

WD6010EM

WD6010M

WD6010T

WD6010TM

WD7010

WD7010M

WD7010T

WD7010TM

Z5000 (art. 900352)

INTERNAL FILTER KIT

- supplementary filter kit, obligatory for dual use machines, operating shoes and surgical instruments



Machines

WD1160	WD1160P4	WD3060	WD4060	WD4060P4
WD4190	WD4190P4	WD5090T	WD5090TFO	

ZCS1 (art. 900042)

UPPER LEVEL TROLLEY FOR OPERATING ROOM SHOE WASHING

- equipped with sprayer arm
- capacity: 4 pairs of standard operating room shoes
- made of stainless steel



Machines

WD1160	WD1160P4	WD3060	WD4060	WD4060P4
WD4190	WD4190P4	WD5090T	WD5090TFO	

ZCS2 (art. 900022)

LOWER LEVEL TROLLEY FOR OPERATING ROOM SHOE WASHING

- level with machine bottom sprayer
- capacity: 8 pairs of standard operating room shoes
- made of stainless steel



Machines

WD1160	WD1160P4	WD3060	WD4060	WD4060P4
WD4190	WD4190P4	WD5090T	WD5090TFO	

ZCS3 (art. 900014)

2-LEVEL WASHING TROLLEY FOR OPERATING ROOM SHOE WASHING

- upper level with sprayer arm
- lower level with machine bottom sprayer
- capacity: 12 pairs of standard operating room shoes
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

Mini-invasive and videolaparoscope trolleys

EF2DS (art. 902155)

RIGID VIDEOLAPAROSCOPE AND ENDOSCOPE TROLLEY WITH DRYING CONNECTOR

- n. 38 injectors of various diameters for endoscope cable instruments (to place when disassembled) es. cannulas, set of rigid endoscopes, videolaparoscope, trocars
- n. 10 special connectors for treating hollow instruments with a max length of 550 mm
- n. 2 baskets for very small silicone instruments and small, different diameter components
- n. 1 basket for placing endoscope mechanical equipment with different diameter silicone supports
- equipped with a connector with drying system coupling
- made of stainless steel



Machines

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

EF3 (art. 901217)

OPHTHALMOLOGY TROLLEY

- n. 12 male connectors diam. 6, 4, 3 mm for hollow instruments and tube accessories
- n. 9 female connectors diam. 7 and 4 mm
- n. 2 spring check valves with respective clamps for external washing of delicate tubular instruments
- n. 5 hollow supports diam. 7 and 4 mm for internal and external washing of cannulas
- n. 4 baskets for very small silicone instruments and small, different diameter components
- n. 2 basket for positioning of tiny instruments with different diameter silicone supports
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

EF3DS (art. 902889)

OPHTHALMOLOGY TROLLEY WITH DRYING CONNECTOR

- n. 12 male connectors diam. 6, 4, 3 mm for hollow instruments and tube accessories
- n. 9 female connectors diam. 7 and 4 mm
- n. 2 spring check valves with respective clamps for external washing of delicate tubular instruments
- n. 5 hollow supports diam. 7 and 4 mm for internal and external washing of cannulas
- n. 4 baskets for very small silicone instruments and small, different diameter components
- n. 2 basket for positioning of tiny instruments with different diameter silicone supports
- made of stainless steel



Machines

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

EF5 (art. 903660)

RIGID VIDEOLAPAROSCOPE AND ENDOSCOPE TROLLEY

- n. 20 injectors with supports for hollow instruments with diam. from 1 to 10 mm
- suitable for treating instruments with a max length of 550 mm
- n. 20 Luer-lock quick connectors
- n. 2 baskets for positioning of small and standard instruments
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

EF5DS (art. 903661)

RIGID VIDEOLAPAROSCOPE AND ENDOSCOPE TROLLEY WITH DRYING CONNECTOR

- n. 20 injectors with supports for hollow instruments with diam. from 1 to 10 mm suitable for treating instruments with a max length of 550 mm
- n. 20 Luer-lock quick connectors
- n. 2 baskets for positioning of small and standard instruments
- equipped with a connector with drying system coupling
- made of stainless steel



Machines

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

Container washing trolleys

CPC1B (art. 903232)

TROLLEY FOR WASHING CONTAINERS

- capacity: 2 containers
- container maximum dimensions LxDxH: 300 x 300 x 150 + 2 covers
- made of stainless steel



Machines

WD1160
WD4190

WD1160P4
WD4190P4

WD3060
WD5090T

WD4060
WD5090TFO

WD4060P4

CPCDIN415 (art. 904372)

TROLLEY FOR WASHING CONTAINERS

- capacity: 4 containers
- container maximum dimensions LxDxH: 600 x 300 x 270 + 4 covers
- n. 22 injection nozzles with detachable tip
- made of stainless steel

CPCENDO (art. 903783)

TROLLEY FOR ENDOSCOPE CONTAINERS

- trolley for washing one endoscope container with cover
- cover maximum dimensions LxDxH: 510 x 410 x 30 mm
- container maximum dimensions LxDxH: 515 x 430 x 120 mm
- to be used together with the ceiling spray arm IRCP
- made of stainless steel



Machines

WD1160
WD4190

WD1160P4
WD4190P4

WD3060
WD5090T

WD4060
WD5090TFO

WD4060P4

SCDIN10 (art. 903406)

CONTAINER WASHING SUPPORT

- capacity: 1 DIN container with lid
- can be placed on rack CSK1210UNIV (max. 2)
- made of stainless steel



Machines

WD6010

WD6010E

WD6010EM

WD6010M

WD6010T

WD6010TM

WD7010

WD7010M

WD7010T

WD7010TM

Anesthesia washing trolleys

AE2BDS (art. 903300)

ANESTHESIA SET WASHING TROLLEY WITH DRYING SYSTEM

- n. 8 injection supports for corrugated tubes max length 1500 mm
- n. 13 nozzles for catheters
- n. 3 balloon nozzles
- n. 4 supports for masks and accessories
- equipped with a connector with drying system coupling (also needed by machines without drying system for right positioning of the instruments)
- made of stainless steel



Machines

WD4060 WD4060P4 WD4190 WD4190P4 WD5090T
WD5090TFO

CSKAE7 (art. 903580)

ANESTHESIA SET WASHING TROLLEY

- n. 8 injection supports for corrugated tubes max length 1500 mm
- n. 10 nozzles for catheters
- n. 4 balloon nozzles
- n. 4 supports for masks and accessories
- made of stainless steel



CSKAE810 (art. 903517)

ANESTHESIA SET WASHING TROLLEY

- n. 8 injection supports for corrugated tubes max length 1500 mm
- n. 10 nozzles for catheters
- n. 4 balloon nozzles
- n. 4 supports for masks and accessories
- made of stainless steel



Machines

WD6010

WD6010E

WD6010EM

WD6010M

WD6010T

WD6010TM

WD7010

WD7010M

WD7010T

WD7010TM

Special trolleys

CM3 (art. 901969)

VARIOUS MASK WASHING TROLLEY

- upper level with sprayer arm
- lower level with machine bottom sprayer
- capacity: 24 standard masks
- made of stainless steel



Machines

WD1160

WD1160P4

WD3060

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

WD5090TFO

CS1CYST (art. 903282)

UROLOGY CYSTOSCOPES RACK

- upper rack with drying connection
- equipped with support for washing 2 cystoscopes (flexible endoscopy) with chemical disinfection cycle
- n. 2 Luer-lock connections each cystoscope for washing channels
- connection for process control (Lum-Check type)
- spray arm for cleaning the external parts



Machines

WD4060

WD4060P4

WD4190

WD4190P4

WD5090T

Transport systems

DIS-WD5 (art. 902411)

RELEASE TROLLEY FOR 60-90 CM MACHINES WITH PLINTH

- for loading and unloading of trolleys to be positioned on the lower level of the machine
- equipped with 4 pivoting and 2 lockable wheels
- equipped with double handle for movement
- lower level for holding supports or other trolleys
- system for collection of liquid residue as required by ISO EN 15883-1 in force
- made of stainless steel



Machines

WD1160	WD1160P4	WD3060	WD4060	WD4060P4
WD4190	WD4190P4	WD5090T	WD5090TFO	

DIS-WD6010 (art. 903455)

RELEASE TROLLEY FOR DOUBLE DOOR MACHINES

- release trolley with safety block system to load/unload washing racks into the washer disinfectors
- diamond-shaped basin for collecting liquid dripping from load material
- equipped with 4 rotating and lockable wheels
- equipped with ergonomic handle for transportation
- safety system for locking washing trolleys
- system for collection of liquid residue as required by ISO EN 15883-1 in force
- made of AISI 304 stainless steel



Machines

WD6010	WD6010E	WD6010EM	WD6010M	WD6010T
WD6010TM				

DIS-WD7 (art. 902218)

RELEASE TROLLEY FOR DOUBLE DOOR MACHINES

- release trolley with safety block system to load/unload washing racks into the washer disinfectant
- diamond-shaped basin for collecting liquid dripping from load material
- equipped with 4 rotating and lockable wheels
- equipped with ergonomic handle for transportation
- safety system for locking washing trolleys
- system for collection of liquid residue as required by ISO EN 15883-1 in force
- made of AISI 304 stainless steel



DIS-WD7010 (art. 903456)

RELEASE TROLLEY FOR DOUBLE DOOR MACHINES

- release trolley with safety block system to load/unload washing racks into the washer disinfectant
- diamond-shaped basin for collecting liquid dripping from load material
- equipped with 4 rotating and lockable wheels
- equipped with ergonomic handle for transportation
- safety system for locking washing trolleys
- system for collection of liquid residue as required by ISO EN 15883-1 in force
- made of AISI 304 stainless steel



Machines

WD7010

WD7010M

WD7010T

WD7010TM

DIS-WD7015 (art. 903862)

RELEASE TROLLEY

- release trolley with safety block system to load/unload washing racks into the washer disinfectant
- diamond-shaped basin for collecting liquid dripping from load material
- equipped with 4 rotating and lockable wheels
- equipped with ergonomic handle for transportation
- safety system for locking washing trolleys
- system for collection of liquid residue as required by ISO EN 15883-1 in force
- made of AISI 304 stainless steel



Machines

WD7015

WD7015T

Supports and baskets

Ultrasound accessories

VU03BIC (art. 902169)

BEAKER FOR ULTRASOUND CLEANER

- volume: 600 ml
- made of PYREX glass



Machines

VU03H

VU03CBIC (art. 902171)

COVER FOR BEAKER 600 ML

- cover made of plastic
- silicone rubber o-ring (cod. VUORING) included
- fit on beaker volume 600 ml



VU03COVER (art. 902099)

3 LITER ULTRASOUND CLEANER COVER

- made of plastic



Machines

VU03H

VU03CSK (art. 902101)

BASKET FOR 3 LITER ULTRASOUND CLEANER

- equipped with lifting handles
- dimensions LxDxH: 198 x 106 x 50 mm
- made of stainless steel



Machines

VU03H

VU03SUP (art. 902165)

BEAKER SUPPORT FOR 3 LITER ULTRASOUND CLEANER

- n. 2 beaker positions
- made of stainless steel



Machines

VU03H

VU06COVER (art. 902157)

6 LITER ULTRASOUND CLEANER COVER

- made of plastic



Machines

VU06H

VU06CSK (art. 902161)

BASKET FOR 6 LITER ULTRASOUND CLEANER

- equipped with lifting handles
- dimensions LxDxH: 255 x 115 x 75 mm
- made of stainless steel



Machines

VU06H

VU06SUP (art. 902166)

BEAKER SUPPORT FOR 6 LITER ULTRASOUND CLEANER

- n. 2 beaker positions
- made of stainless steel



Machines

VU06H

VU13COVER (art. 902100)

13 LITER ULTRASOUND CLEANER COVER

- made of plastic



Machines

VU13H

VU13CSK (art. 902102)

BASKET FOR 13 LITER ULTRASOUND CLEANER

- equipped with lifting handles
- dimensions LxDxH: 250 x 190 x 115 mm
- made of stainless steel



Machines

VU13H

VU13SUP (art. 902167)

BEAKER SUPPORT FOR 13 LITER ULTRASOUND CLEANER

- n. 2 beaker positions
- made of stainless steel



Machines

VU13H

VU30COVER (art. 902158)

30 LITERS ULTRASOUND CLEANER COVER

- made of plastic



Machines

VU30H

VU30CSK (art. 902162)

BASKET FOR 30 LITER ULTRASOUND CLEANER

- equipped with lifting handles
- dimensions LxDxH: 455 x 250 x 115 mm
- made of stainless steel



Machines

VU30H

VU30SUP (art. 902168)

BEAKER SUPPORT FOR 30 LITER ULTRASOUND CLEANER

- n. 2 beaker positions
- made of stainless steel



Machines

VU30H

VU45COVER (art. 902159)

45 LITER ULTRASOUND CLEANER COVER

- made of stainless steel



Machines

VU45H

VU45CSK (art. 902163)

BASKET FOR 45 LITER ULTRASOUND CLEANER

- equipped with lifting handles
- dimensions LxDxH: 455 x 270 x 194 mm
- made of stainless steel



Machines

VU45H

VU90COVER (art. 902160)

90 LITER ULTRASOUND CLEANER COVER

- made of stainless steel



Machines

VU90H

VU90CSK (art. 902164)

BASKET FOR 90 LITER ULTRASOUND CLEANER

- equipped with lifting handles
- dimensions LxDxH: 545 x 450 x 250 mm
- made of stainless steel



Machines

VU90H

VUMICRO1 (art. 902470)

**TINY SURGICAL INSTRUMENT BASKET IN MICRO-STRETCHED
FOR ULTRASOUND CLEANING - DIAMETER 59MM**

- suitable for holding tiny instruments and small instruments in general
- dimensions DIAMETER: 59 mm
- mesh 2 x 2 mm
- made of stainless steel



Thermal sealer accessories

SGROL (art. 902225)

CUTTING DEVICE FOR STERILIZATION POUCH ROLLS

- cutting length 750 mm
- dimensions LxDxH: 873 x 420 x 125 mm
- made of stainless steel



SGROL-DSP (art. 902226)

CUTTING DEVICE FOR STERILIZATION POUCH ROLLS WITH ROLL HOLDER

- pouch film rolls max. 300 mm
- cutting length 750 mm
- dimensions LxDxH: 873 x 425 x 425 mm
- made of stainless steel



SGTRAR (art. 902085)

ROLLER FOR STERILIZATION POUCH SLIDING

- for rotating thermal sealer series SGR only
- dimensions LxDxH: 830 x 280 x 60 mm



SGTRAY (art. 902086)

TRAY FOR STERILIZATION POUCH SLIDING

- for rotating thermal sealer series SGR only
- dimensions LxDxH: 850 x 280 x 60 mm
- made of stainless steel

Special supports

BFS (art. 902881)

SPECULUM SUPPORT

- support for n. 6 single speculum
- dimensions LxDxH: 170 x 275 x 160 mm
- made of stainless steel



CSKFLEX (art. 902884)

WASHING BASKETS FOR HOLLOW INSTRUMENTS

- composed by 2 DIN baskets with direct injection system
- capacity: 10 hollow instruments inserts 400 mm length
- Luer-lock connectors included
- made of AISI 304 stainless steel



INSPMIC10 (art. 904273)

MINI-INVASIVE SURGERY TROLLEY

- n. 5 ULP1
- n. 5 ULBEF2
- n. 5 Luer-lock connections
- n. 5 MIC connection
- n. 8 vertical nozzles
- suitable for trolleys CSKxx10
- made in stainless steel

INSPMIC15 (art. 904274)

MINI-INVASIVE SURGERY TROLLEY

- n. 7 ULP1
- n. 7 ULBEF2
- n. 6 Luer-lock connections
- n. 6 MIC conenction
- n. 10 vertical nozzles
- suitable for trolleys CSKxx15
- realizzato in acciaio inox

KITMIC10 (art. 903499)

CANNULATED INSTRUMENTS BASKET

- CSKDIN basket suitable for MIC / MIS instruments - injection wash
- n. 4 Luer-lock connectors
- n. 5 hole connectors



Machines

WD6010

WD6010E

WD6010EM

WD6010M

WD6010T

WD6010TM

WD7010

WD7010M

WD7010T

WD7010TM

Baskets

CB18 (art. 900016)

BABY BOTTLE CONTAINER FOR 18 250 ML BOTTLES

- dimensions LxDxH: 210 x 425 x 195 mm
- made of stainless steel



CB18/100 (art. 900273)

BABY BOTTLE CONTAINER FOR 18 100 ML BOTTLES

- dimensions LxDxH: 220 x 420 x 115 mm
- made of stainless steel



CB18/120 (art. 900341)

BABY BOTTLE CONTAINER FOR 18 120 ML BOTTLES

- dimensions LxDxH: 220 x 420 x 125 mm
- made of stainless steel



CPS20 (art. 900824)

BASKET WITH LID

- dimensions LxDxH: 430 x 190 x 60 mm
- made of stainless steel



CSK1 (art. 900019)

INSTRUMENT BASKET WITH HANDLES

- suitable for holding of large-sized instruments
- dimensions LxDxH: 450 x 450 x 55 mm
- mesh 5 x 5 mm
- made of stainless steel



CSK1/3 (art. 900522)

INSTRUMENT BASKET WITH HANDLES

- suitable for holding various medium-sized instruments
- dimensions LxDxH: 435 x 160 x 50 mm
- mesh 5 x 5 mm
- made of stainless steel



CSK1/6 (art. 901264)

INSTRUMENT BASKET WITH COVER

- suitable for holding various small-dimension instruments
- dimensions LxDxH: 130 x 225 x 50 mm
- cover mesh 5 x 5 mm
- basket mesh 3 x 3 mm
- made of stainless steel



CSK2 (art. 900020)

BASKET WITH HANDLES


- suitable for holding medium-sized instruments
- dimensions LxDxH: 450 x 225 x 50 mm
- mesh 5 x 5 mm
- made of stainless steel



CSKDIN (art. 900018)

INSTRUMENT BASKET WITH HANDLES

- suitable for holding medium-sized instruments
- dimensions LxDxH: 480 x 245 x 50 mm
- mesh 5 x 5 mm
- made of stainless steel

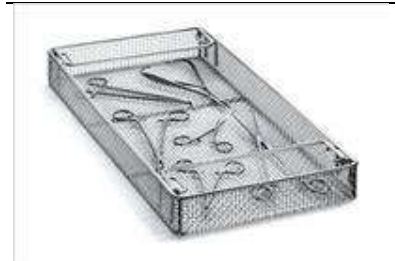

CSKL (art. 900751)

BASKET FOR SURGICAL INSTRUMENTS WITH HANDLES

- dimensions LxDxH: 280 x 600 x 55 mm
- space needed half level of CBS1 and CBS2
- mesh 5 x 5 mm
- made of stainless steel



CSMICRO1 (art. 900811)

MICRO-INSTRUMENTS BASKET IN FINE-DRAWN MESH WITH COVER

- suitable for holding tiny and small instruments
- dimensions LxDxH: 240 x 235 x 55 mm
- 2 x 2 mm mesh
- made of stainless steel



CSMICRO2 (art. 901234)

MICRO-INSTRUMENTS BASKET IN FINE-DRAWN MESH WITH COVER

- suitable for holding tiny and small instruments
- dimensions LxDxH: 100 x 80 x 45 mm
- mesh 2 x 2 mm
- made of stainless steel



CSMICRO3 (art. 902056)

MICRO-INSTRUMENTS BASKET IN FINE-DRAWN MESH WITH COVER

- suitable for holding tiny and small instruments in general
- dimensions LxDxH: 280 x 168 x 45 mm
- mesh 2 x 2 mm
- made of stainless steel



D-CM1 (art. 901265)

TALL BASKET WITH HANDLES

- suitable for holding various instruments
- dimensions LxDxH: 300 x 155 x 115 mm
- mesh 3 x 3 mm
- made of stainless steel



D-SV (art. 901200)

SUPPORT FOR VERTICAL INSTRUMENTS

- suitable for holding various vertical instruments
- dimensions LxDxH: 300 x 155 x 150 mm
- made of stainless steel



SA3 (art. 900536)

BASIN AND KIDNEY WASHING SUPPORT


- suitable for holding kidney bowls and basins
- dimensions LxDxH: 265 x 170 x 65 mm
- made of stainless steel



SBAS5 (art. 904183)

SUPPORT FOR 5 BOWLS


- dimensions LxDxH: 135 x 490 x 60 mm
- support for positioning of n. 5 bowls (max diameter: 345 mm, max depth: 120 mm)

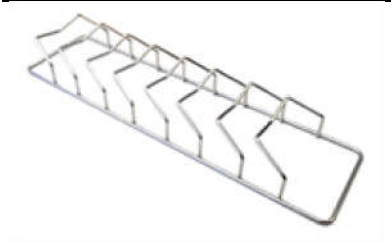

SBAS6 (art. 904125)

SUPPORT FOR 6 BOWLS

- dimensions LxDxH: 475 x 485 x 170 mm
- support for positioning of n. 6 bowls (max diameter: 320 mm, max depth: 88 mm)

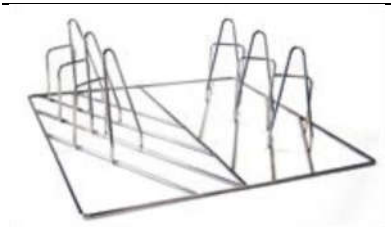

SCVO (art. 903865)

SUPPORT FOR OPTICAL VIDEO CABLES

- dimensioni LxPxH: 110 x 110 x 200 mm



Water systems accessories

PD50 (art. 900984)

AUXILIARY TANK 50L

- 50 liters tank
- suitable to provide demineralized water to the washer machine
- made of PE-HD
- 3/4" BSP tap
- dimensions DIAMETERxH: 340 x 640 mm



SDA100 (art. 901021)

AUXILIARY TANK 100L WITH AUTOMATIC LEVEL CONTROL

- 100 liters tank
- suitable to provide demineralized water to the washer machine
- made of PE-HD
- with automatic level control for external water feeding
- compatible with demi water booster pump PAD / PAD2
- rapid connection for 8 mm pipes in PE
- connection for demineralised water inlet
- connection for water inlet by auxiliary system
- 3/4" BSP tap
- dimensions DIAMETERxH: 450 x 820 mm



SDA50 (art. 901020)

AUXILIARY TANK 50L WITH AUTOMATIC LEVEL CONTROL

- 50 liters tank
- suitable to provide demineralized water to the washer machine
- made of PE-HD
- with automatic level control for external water feeding
- compatible with demi water booster pump PAD / PAD2
- rapid connection for 8 mm pipes in PE
- connection for demineralised water inlet
- connection for water inlet by auxiliary system
- 3/4" BSP tap
- dimensions DIAMETERxH: 340 x 640 mm



Consumables

CARTADP24 (art. 901367)

PRINTER PAPER DP24/3RTC

- 10 roll pack



INKDP24 (art. 901368)

INK FOR PRINTER DP24/3RTC

- single package cartridge

KTANK (art. 904162)

TOOL FOR OPENING JERRYCANS

- double-key for opening Dr. Weigert and SMEG canisters



SACCRES (art. 901378)

RESINS DISPOSAL BAG

- exhausted resin disposal bag

SGRINK (art. 902192)

INK FOR THERMAL SEALER PRINTER

- single package cartridge



WD-PAPER (art. 901709)

THERMAL PAPER FOR PRINTER

- package of 10 rolls
- to be used with printers WD-PRINT, WD-PRINTE



WD-RIBBON (art. 904008)

RIBBON FOR WD-PRINTEINK

- ribbon for WD-PRINTEINK
- dimensions 57 mm x 74 M
- n. 4 pieces per box



CL3000 (art. 900186)

EMPTY RESINS CONTAINER

- capacity: 15 liters for WP3000
- made of stainless steel AISI 316L
- jolly fast connectors included



MI700WP (art. 900166)

RESIN KIT REFILL - FOR 15 LITER CONTAINER

- mixed bed resins kit for WP3000 water purifier
- supplied with used resin disposal bag
- stainless steel resin container CL3000 not included



980103 (art. 980103)

SPECIAL SALT TABLETS FOR EXTERNAL WATER SOFTENERS

- 25 Kg pack



SALE (art. 901358)

SPECIAL SALT FOR INCORPORATED WATER SOFTENER

- single pack of 20 kg (20 boxes of 1 Kg each)

SGCHECK (art. 902089)

SEALING QUALITY TEST

- tester for sealing quality verification
- pack of 250 testers



Detergents and additives

ACIDGLC2 (art. 901117)

ACIDGLASS C2

- liquid phosphate-free neutralizer
- containing citric acid monohydrate
- 5 liters can



DETERLQC2 (art. 901116)

DETERLIQUID C2

- strong universal liquid alkaline detergent phosphate-free for washing steel up to 95 °C
- also suitable for washing rubber and plastic at 60 °C
- containing sodium silicate, polyacrylate, potassium hydroxide
- 5 liters can



ECOSEPT (art. 903626)

ECOSEPT

- liquid disinfection additive for disinfection of delicate materials (ex. operating room shoes)
- containing 2-PHENOXYETHANOL
- 5 liter can

