

# Haier Biomedical

Intelligent Protection of Life Science

## -86°C ULT Freezer

DW-86L490J



### Energy saving

VIP(Vacuum Insulation Panel) design with extra insulation cuts down on heat loss to cabinet by 25% and helps reduce energy consumption by more than 20%, power consumption only 11.5 kWh/day



### Double door design

More convenient to be used by different users



### Reliable

Cabinet temperature can reach -86°C with temp uniformity within 4°C

### Complementary products



CO<sub>2</sub>  
Backup  
Systems



Chart  
Recorder



Stainless  
Steel  
Racks



Cryo  
Boxes

# Specifications

## DW-86L490J

### Technical Data

Type	-86°C ULT Freezer
Cabinet Type	Upright
Climate Class	N
Cooling Type	Direct cooling
Defrost Mode	Manual
Refrigerant	HC
Noise Level (dB)	50

### Performance

Cooling performance(°C)	-86
Temperature Range(°C)	-40~-86

### Control

Controller	Microprocessor
Display	LED

### Electrical Data

Power Supply(V/Hz)	220-240/50	208-230/60
Power(W)	900	1000
Electrical Current(A)	8	8

### Functions

High/Low Temperature	Y
Remote Alarm	Y
Hot Condenser	Y
Power Failure	Y
Sensor Error	Y
Low Battery	Y
High Ambient Temp	Y
Door Ajar	Y



### Accessories

Caster	Y
Foot	Y
Test Hole	Y/2
Shelves/Inner Doors	3/4
Temperature Recorder	Optional
USB Interface	Y
Rs232/485 Interface	Optional
CO <sub>2</sub> Backup System	Optional
LN <sub>2</sub> Backup System	Optional
Racks & Boxes	Optional

### Dimensions

Capacity(L/Cu.Ft)	490/17.3	
Net/Gross Weight(approx)	295/335(kg)	650.4/738.5(lbs)
Interior Dimensions(W*D*H)	590*630*1310(mm)	23.2*24.8*51.6(in)
Exterior Dimensions(W*D*H)	860*900*1980(mm)	33.9*35.4*78.0(in)
Packing Dimensions(W*D*H)	925*985*2150(mm)	36.4*38.8*84.6(in)
Container load (20'/40'/40'H)	12/24/24	

# Accessories & consumables

## DW-86L490J

### Consumables

Racks	Part Code	
Model	DCJ-54-A	DCJ-54-B (the 1 <sup>st</sup> and the 2 <sup>nd</sup> layer)
	DCJ-44-A	DCJ-44-B (the 3 <sup>rd</sup> and the 4 <sup>th</sup> layer)
Q'ty/Layer	4	
Dimension (W*D*H)	138*565*293(mm)	140*565*293(mm)
	5.4*22.2*11.5(in)	5.5*22.2*11.5(in)
Dimension (W*D*H)	138*565*240(mm)	140*565*240(mm)
	5.4*22.2*9.4(in)	5.5*22.2*9.4(in)

Boxes	Part Code			
Model	T130	HBB-200-PP	HBB-281-PC	HBB-200-C
Material	PP	PP	PC	PC
Dimension (W*D*H)	133*133*53.5(mm) 5.2*5.2*2.1(in)			
Q'ty/Rack	20			
	16			
<b>Total Q'ty</b>				
Racks	8			
Boxes	160			
	128			
Vials	12960	16000	12960	16000
	10368	12800	10368	12800

### Complementary products



CO2 Backup Systems



Stainless Steel Racks


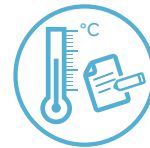


Chart Recorder

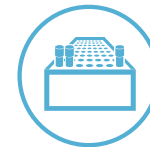

Cryo Boxes

### CO<sub>2</sub> Backup System

Model	HBX-IC
Temperature range	-40°C to -70°C

### LN<sub>2</sub> Backup System

Model	HBX- II A
Temperature range	-90°C to -135°C

### Chart Recorder

Temperature range	+50°C to -115°C
-------------------	-----------------

### Others

Certificate	CE /
-------------	------



Scan 2D barcode with your mobile device to learn more information about Haier Biomedical

#### Haier Biomedical UK

Address: 2 Birch House, Brambleside, Bellbrook Industrial Estate, Uckfield, East Sussex, TN22 1QQ, United Kingdom  
 TEL: +44(0) 1825 744 696 (Head Office)  
[www.haierbiomedical.co.uk](http://www.haierbiomedical.co.uk)

Intelligent Protection  
of Life Science